

Samson Power Drive PD-6A[®]



Advanced Technologies That Move You

The adjustable PD-6A handlebars are ergonomically designed for comfort and easy operation. The PD-6A exclusive dynamic brake system takes away the fear and physical stress of steep slopes.



CAR



TRAIN




SHIP



AIRPLANE



STAIRS

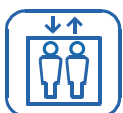

ELEVATOR



TZORA
Active Systems Ltd

Samson Power Drive PD-6A®

INFUSE YOUR WHEELCHAIR WITH ENERGY!

The battery powered PD 6A attaches to any manual wheelchair in seconds and takes on the chore of pushing the chair. Powerful enough for steep inclines and long range hauls, the PD-6A is extremely lightweight, easy to handle, easy to attach and detach. It's ideal for transportation in the trunk of a car.

- Easily clamps on and detaches from wheelchair.
- Forward and backward operation.
- Long range and plenty of power with the 12 A/h battery.
- Adjustable height, ergonomically designed handlebars.
- Dynamic, electromagnetic brake system.
- Battery is easily removed for recharging.
- On/off button on left and right handles.
- Forward and reverse control on panel.

Technical Specifications:

Carrying Capacity	115 kg	253 lbs
Weight:		
Total Weight (12A/h)	19 kg	42 lbs
Weight of PD-6 without batteries	10 kg	22 lbs
Weight of the battery 12A/h	9 kg	20 lbs
Charger	100-240	VAC
Tiller Type Adgustable		
Speed Adjustable up to	6 km/h	4 mile/h
Range		
12 Ah battery	up to 15 km	up to 9 mile

Weights & measurements are approximate



PD-6A



PD-6A
Control
Panel



PD-6A
attached to
a wheelchair

